

Venturi Barrels for Intersurgical Masks

24% (Blue) Venturi Barrel – Intersurgical

- Recommended start flow rate = 2 L/MIN

28% (White) Venturi Barrel – Intersurgical

- Recommended start flow rate = 4 L/MIN

31% (Orange) Venturi Barrel - Intersurgical

- Recommended start flow rate = 6 L/MIN

35% (Yellow) Venturi Barrel – Intersurgical


- Recommended start flow rate = 8 L/MIN

40% (Red) Venturi Barrel – Intersurgical


- Recommended start flow rate = 10 L/MIN

60% (Green) Venturi Barrel – Intersurgical

- Recommended start flow rate = 15 L/MIN

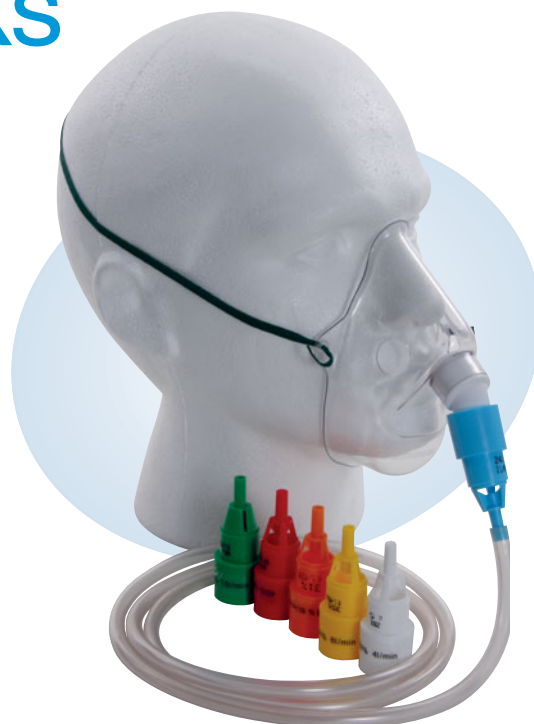

Venturi Barrels for Intersurgical Masks

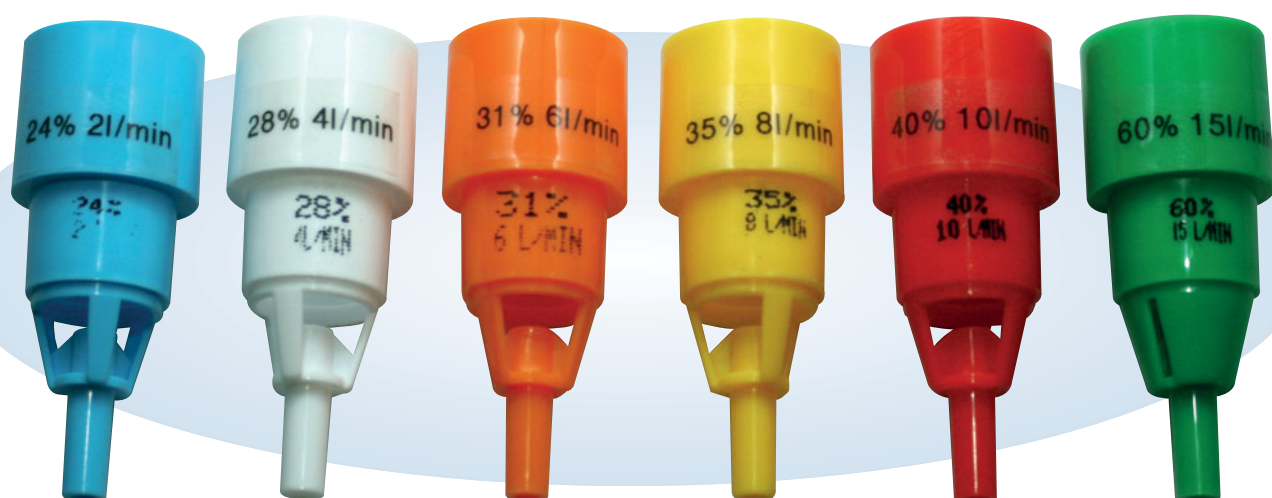

Venturi Barrels for Intersurgical Masks